


ROB Knee Holder System

The **ROB Knee Holder System** is especially designed for all standard TKA applications and robotic or navigated training procedures. The femoral head support & hemisphere set is mimicking the hip joint for a full anatomical range of motion. This system holds various SYNBONE® Knee models with notch in a surgical position. This compact system guarantees a standardized education concept and convinces in quality and performance. The unique notch holders prevent the bones from slipping or moving while surgeons practice various procedures for the total knee replacement with robotic or navigated training procedures.

Features and benefits

- For robotic or navigated training procedures
- Designed for left and right sided models
- One adapter supports either tibia or femur
- No slipping or moving of the bones
- Firm hold during the exercise
- Tibial external alignment thanks to soft tissue foot
- Maintenance free / easy to clean
- Delivered in a solid transportation case

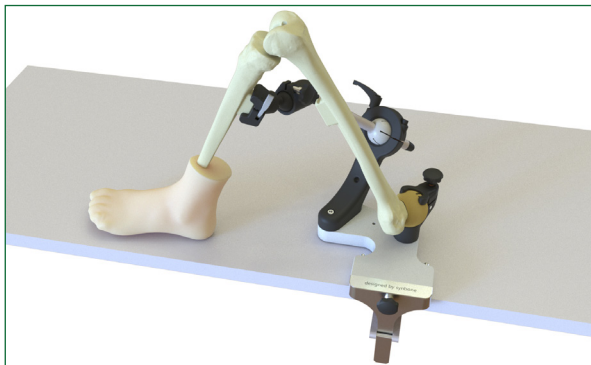


ROB Knee Holder System

Description	ROB Knee Holder System for Total Knee Arthroplasty workshops		
Dimensions	Length:	505 mm	
	Width:	297 mm	
	Height:	610 mm	
Material	Anodised aluminium, stainless steel		
Weight	3.2 kg		



ROB Knee Holder System Case

Description	Solid and waterproof transportation case with foam inlays		
Dimensions	Length:	525 mm	
	Width:	450 mm	
	Height:	235 mm	
Material	PPC Polypropylene Copolymer		
Weight	5.4 kg		

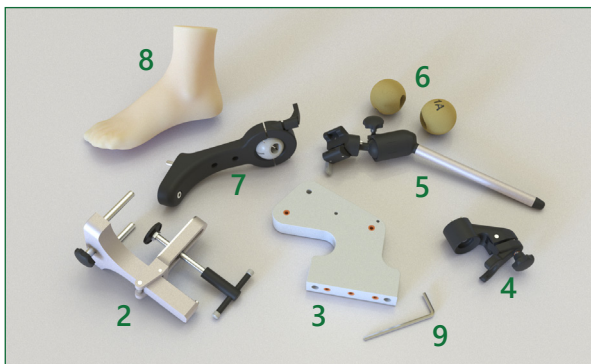


Scope of Delivery


- 1 ROB Knee Holder System Case (optional)
- 2 Table Clamp
- 3 Base Plate
- 4 Femoral Head Support
- 5 Guiding Rail with Knee Support Notch Holder
- 6 Femoral Head Hemisphere Set (2 pcs)
- 7 Knee Support Base
- 8 Reference Foot
- 9 Hex wrench

PRODUCT DETAILS

SYNBONE®
better education outcome

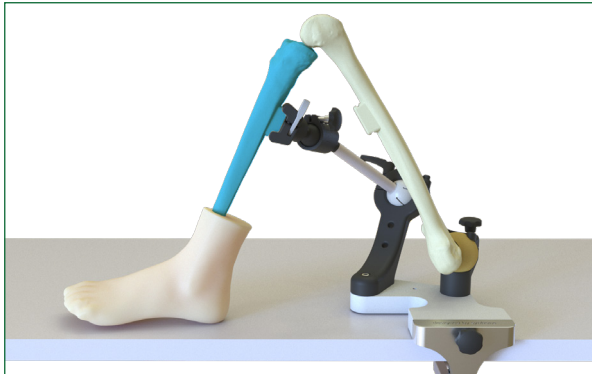


Notch holder securing the tibia (for procedures that focus on the tibia)

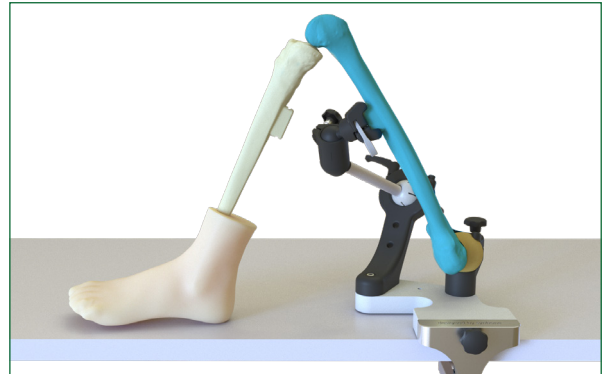


Notch holder securing the femur (for procedures that focus on the femur)



Range of Motion

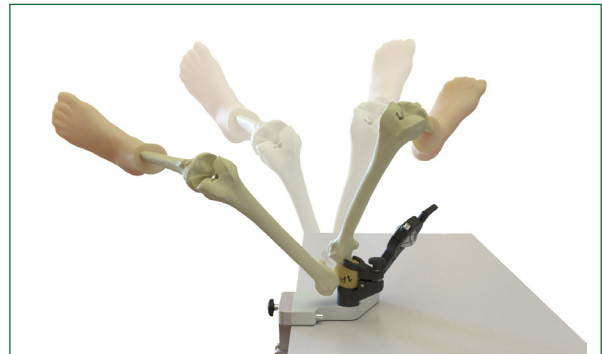

Range of Motion

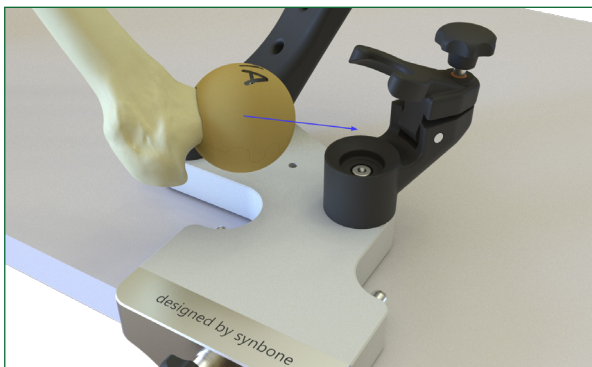

Femur with mounted hemisphere set



Secure attachment to the femoral head support

To obtain best training results we recommend using our SYNBONE® Knee with notch. Please contact us for our complete product portfolio for your Total Knee Arthroplasty trainings.

SYNBONE® is your preferred partner for better education outcome.

Worldwide contact information – sales@synbone.com

contact your
Regional Sales Representative
SYNBONE.COM