


Skull Holder w/ Ball Head

The SYNBONE® Skull Holder w/ Ball Head is the small brother of the SYNBONE® Skull Holder System and is recommended to be used in workshop rooms with a high number of participants. The **Skull Holder with Ball Head** is specially designed for CMF Trauma, Orthognatic and Neuro Education and holds SYNBONE® skulls, partial skulls and maxillas with or without vice attachment. The included 3-pin adapter offers a wide range of applications and a 360° stepless rotation of the skulls.


Features and benefits

- Easy and fast setup thanks to table clamp
- Includes 3-pin adapter with ball head
- 360 ° stepless rotation of models
- Firm hold on working table during the exercise
- No slipping or moving of the skull models
- Maintenance free / easy to clean

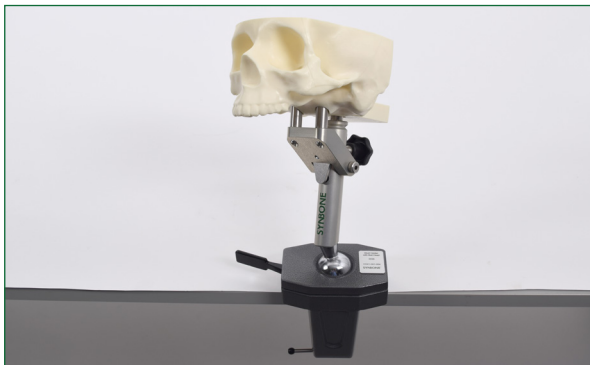

Skull Holder w/ Ball Head

Description	Skull Holder with Ball Head for CMF Trauma, Orthognatic and Neuro workshops
Dimensions	Length: 390 mm Width: 175 mm Height: 80 mm
Material	Anodised aluminium, stainless steel
Weight	1.450 kg

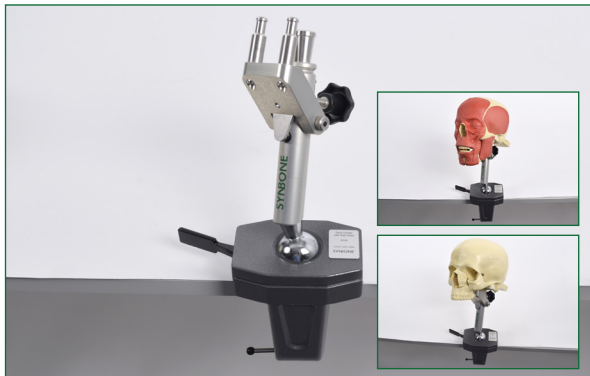

Fast set-up


The **Skull Holder with Ball Head** can be easily mounted on the working table and is ready in seconds. The base of the supporting column is equipped with a stainless steel ball head that allows a 360° stepless adjustment of the model in all directions.

3-pin adapter

The 3-pin adapter holds SYNBONE® skulls, partial skulls and maxillas **with or without vice attachment** securely in place. The three pins are designed to fit in all models and guarantee a firm hold.

Note:

For all Hemi skull applications, please use SYNBONE® Skull Holder System Art. 0038.50 with the additional Hemi skull adapter.

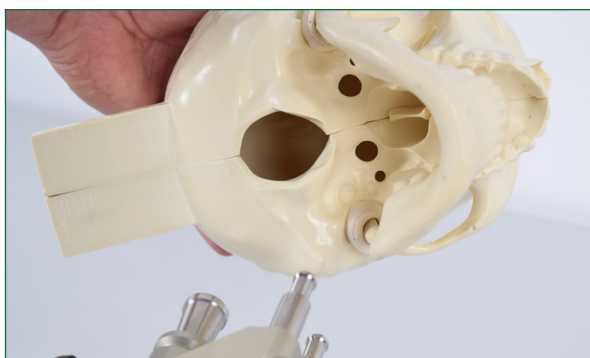

contact
sales@synbone.com